


Flexi-Q EAI

2-Step Disposable Auto-Injector
for emergency use

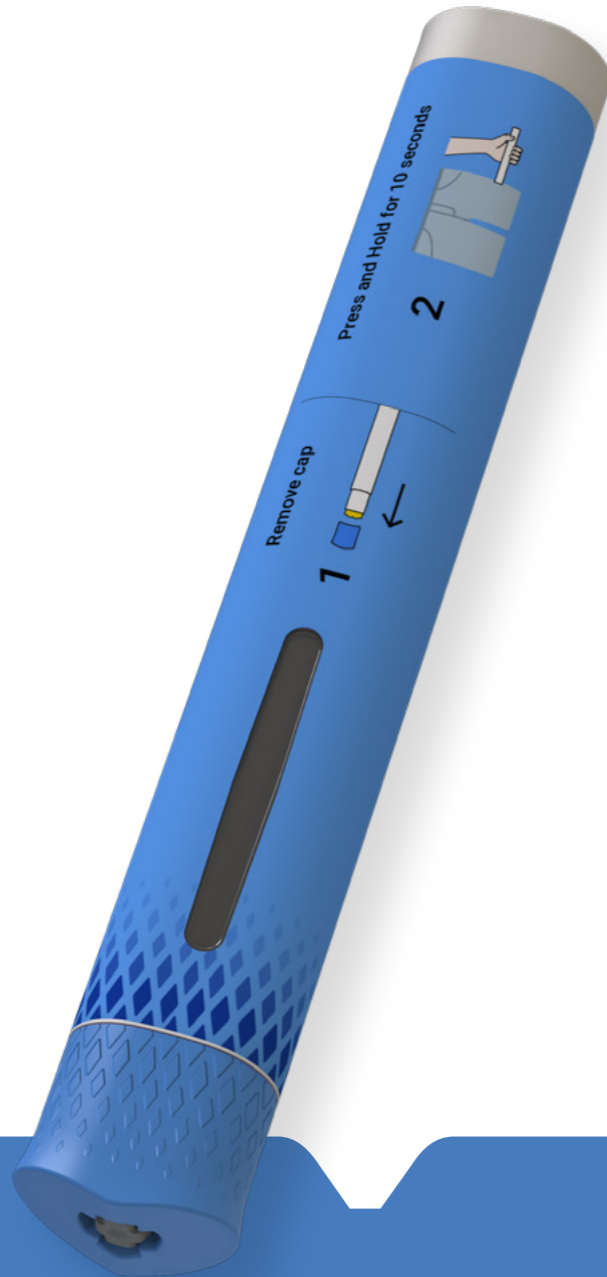

The Flexi-Q EAI is a 2-step fully disposable auto-injector for drugs in 1 mL or 2.25 mL pre-filled syringes (PFS). This easy to use and single-use auto-injector features passive needle shielding and allows better patient safety.

Pharmaceutical Company Advantages

Differentiation

- Flexible delivery volumes, fixed dose:
0.1 - 1.0 mL for the 1 mL Flexi-Q EAI 100 (for 1 mL PFS)
0.1 - 2.25 mL for the 2.25 mL Flexi-Q EAI 225 (for 2.25 mL PFS)
- Compatible with Intramuscular (IM) and Subcutaneous (SC) injections
- Proprietary damping mechanism - preventing impact on glass syringe -
- allowing use of extremely strong injection springs:
a) for high-viscosity drugs
b) for fast injections required for emergency drugs
- Reduced skin reactions - injection starts only after full needle penetration
- Delivery of the entire dose - same dead space as with manual techniques
- Optional large observation window - allows view of the entire syringe fill-volume and injection progress
- Easy and flexible customization options

Short regulatory path

- Compatible with standard glass syringes - no need to change drug container


Patient Benefits

Ease of use

- 2-step activation with fully automatic injection process
- Simple and intuitive operation

Patient safety

- Prevents needle-stick injuries - Passive Needle safety throughout the injection process
- Reduced local skin reactions - injects only into the correct depth (IM / SC)
- Clear visual and audible confirmation of the injection progress from activation to completion
- Fast and simple activation design: -Activation by pressing against the skin - e.g. for emergency use

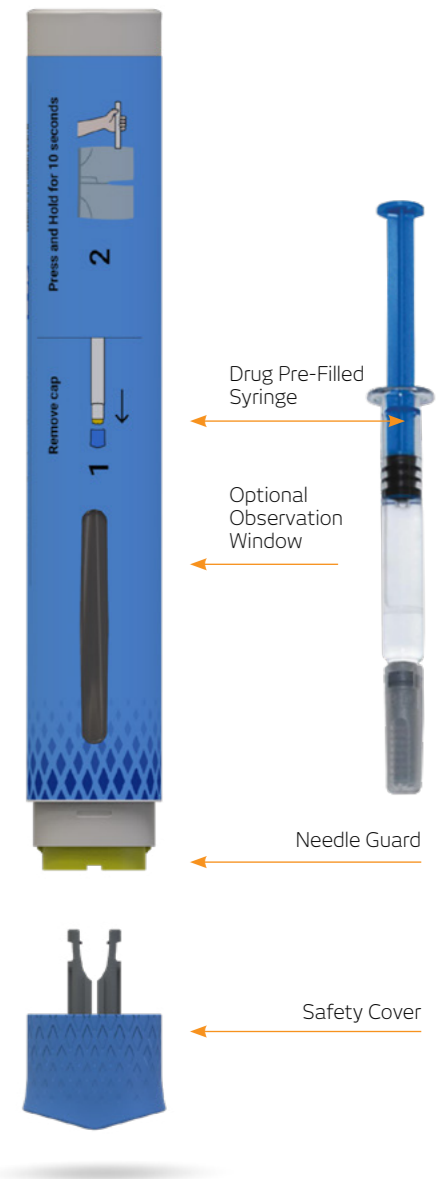

See how easy it is to use:



1. Remove safety cap



2. Push against injection site

* When the auto-injector is removed from the injection site, a Needle Guard automatically covers and locks over the needle for safe disposal of the auto-injector.